


Upload and Manage files to be analyzed with CIBERSORTx

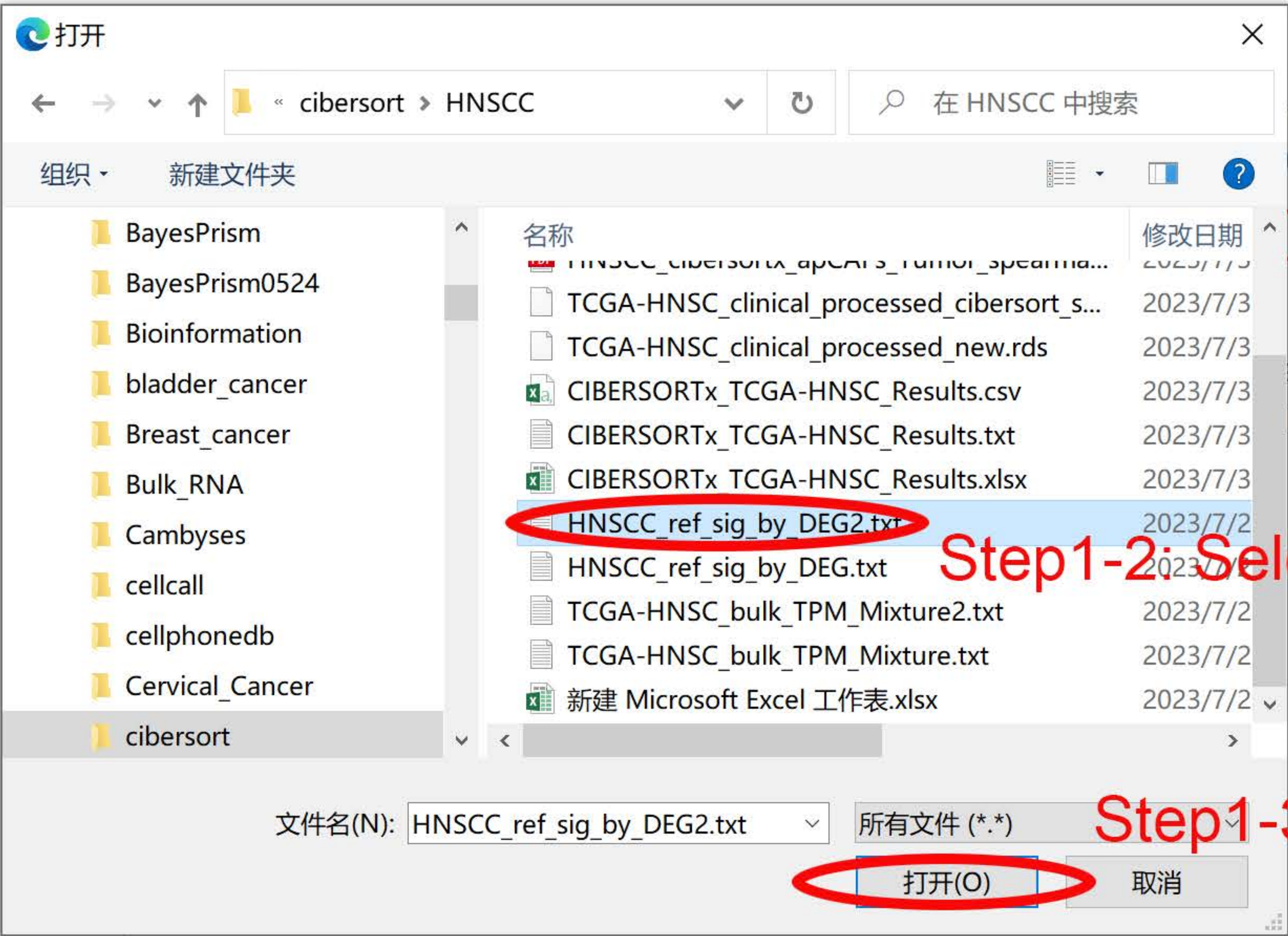

Step1-2: Select signature matrix.

Step1-3: Confirm signature matrix.

Add files...

Start upload

Cancel upload

Delete

☐ Select all

Step1-1: Add signature matrix.

Title	File Type	File name	Size	Upload Date	Actions
BRCA_p2_Mixture	Mixture	TCGA-BRCA_bulk_TPM_p2_Mixture2.txt	70.49 MB	07/04/23 08:12:02 AM	<div>Delete</div>

Instructions:

- Be sure all input files to be uploaded are text tab-delimited and do not contain any double quotation marks. File extension should end in **txt** or **tsv**, otherwise the file will be rejected.
Please note that CSV files are NOT supported; please convert them to text tab-delimited format.
- **File names with white spaces and any special characters (e.g. ~ ` ! @ # \$ % ^ & * () { | } / \ : ; " ' < > ? , = +) are currently not supported. Dashes and underscores (- _) are supported. Please ensure there are no white spaces in the file names.**
- Click 'Add files' or drag and drop files into the browser window.
- Enter a Title for each file, which will be used as an alias for the uploaded file. If no title is entered, the file name will be used.
- Also specify a file type for each file. This will determine which file field the file will appear in the CIBERSORTx form. If no field is specified, the file will **not** appear in the form.
- Uploaded files appear below in reverse chronological order and can be individually deleted
- Note: there is a file quota of 1000 MB. If your file exceeds your available quota, the file will not be saved to the server. (1 MB = 1000000 bytes)
- Please refresh the page after uploads to update quota status.

Quota: 618.0 of 1000 MB used, 382.0 MB remaining

File Type Reference (more...)

Add files...

Start upload

Cancel upload

Delete


☐ Select all

Title	File name	Size	Upload Date	Actions
<div><div>File Type</div><div><div>Title:</div><div>HNSCC_Signature</div><div>File Type:</div><div>Signature Matrix</div></div></div>	HNSCC_ref_sig_by_DEG2.txt	160.90 KB		<div><div>Start</div><div>Cancel</div></div>

Step2-1: Annotate the signature matrix.

Step2-2: Upload the signature matrix.

Upload and Manage files to be analyzed with CIBERSORTx

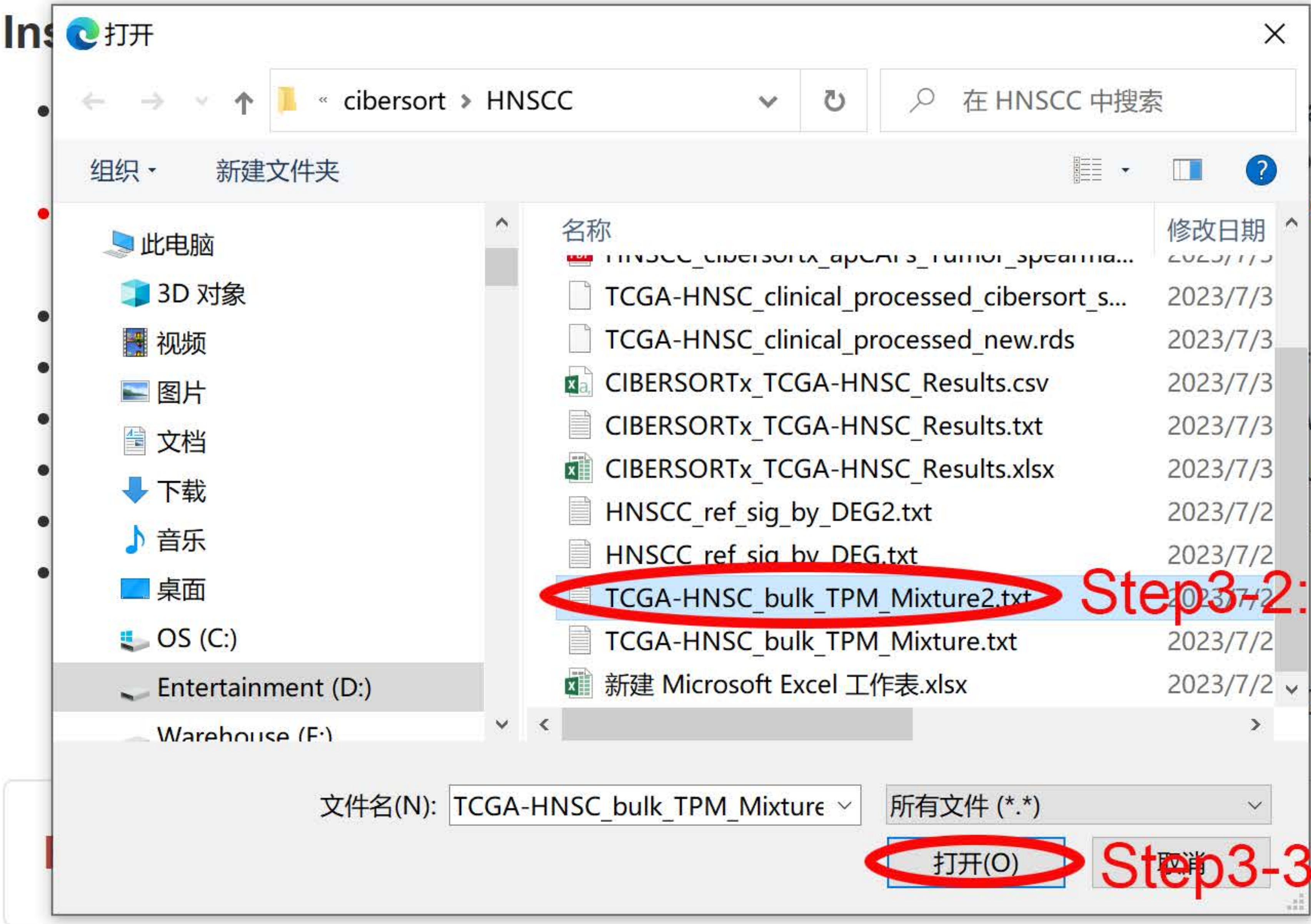

Step3-2: Select Mixture.

Step3-3: Confirm Mixture.

Step3-1: Add Mixture.

+

Add files...

i

Start upload

o

Cancel upload

🗑


Delete

☐ Select all

Title	File name	Size	Upload Date	Actions
HNSCC_Signature	HNSCC_ref_sig_by_DEG2_1.txt	160.90 KB	07/18/23 12:56:56 PM	<div><div>🗑 Delete</div><div><input type="checkbox"/></div></div>
Signature Genes				

Instructions:

- Be sure all input files to be uploaded are text tab-delimited and do not contain any double quotation marks. File extension should end in **txt** or **tsv**, otherwise the file will be rejected.
Please note that CSV files are NOT supported; please convert them to text tab-delimited format.
- **File names with white spaces and any special characters (e.g. ~ ` ! @ # \$ % ^ & * () { | } / \ : ; ' " < > ? , = +) are currently not supported. Dashes and underscores (- _) are supported. Please ensure there are no white spaces in the file names.**
- Click 'Add files' or drag and drop files into the browser window.
- Enter a Title for each file, which will be used as an alias for the uploaded file. If no title is entered, the file name will be used.
- Also specify a file type for each file. This will determine which file field the file will appear in the CIBERSORTx form. If no field is specified, the file will **not** appear in the form.
- Uploaded files appear below in reverse chronological order and can be individually deleted
- Note: there is a file quota of 1000 MB. If your file exceeds your available quota, the file will not be saved to the server. (1 MB = 1000000 bytes)
- Please refresh the page after uploads to update quota status.

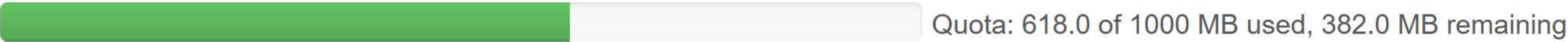

File Type Reference (more...)

+

Add files...

i

Start upload

o

Cancel upload

🗑

Delete

☐ Select all

Title	File name	Size	Upload Date	Actions
<div><div>File Type</div><div><div>Title:</div><div>HNSCC_Mixture</div><div>File Type:</div><div>Mixture</div></div></div>	TCGA-HNSC_bulk_TPM_Mixture2.txt	63.12 MB	<div><div><div>i</div>Start</div><div><div>o</div>Cancel</div></div>	
HNSCC_Signature	HNSCC_ref_sig_by_DEG2_1.txt	160.90 KB	07/18/23 12:56:56 PM	<div><div>🗑</div>Delete</div> <div><input type="checkbox"/></div>

Step4-2: Upload the Mixture.

Step4-1: Annotate the Mixture.

CIBERSORTx: Digital Cytometry



Select Analysis Module: Step5-1: Select impute cell fractions module.

1. Create Signature Matrix

2. Impute Cell Fractions

3. Impute Cell Expression

Currently selected module: 2. Impute Cell Fractions

Select Analysis Mode:

Example

Custom

Step5-2: Select custom mode.

Configure Custom Input:

Signature matrix file*

HNSCC_Signature

Step5-3: Select the uploaded signature matrix.

[More Information...](#)

Mixture file*

HNSCC_Mixture

Step5-4: Select the uploaded Mixture.

[More Information...](#)

☐ Enable batch correction

☒ Disable quantile normalization (disabling is recommended for [RNA-Seq data](#))

☐ Run in absolute mode

Step5-5: Select this option for RNA-seq data.

[More Information...](#)

Permutations for
significance analysis

none

Step5-6: Default parameters.

[More Information...](#)

Run

Step5-7: Run CIBERSORTx.

CIBERSORTx: Digital Cytometry

CIBERSORTx

Run Again

CIBERSORTx Runtime Console

Step6: Download the results.

Download:

TXT

CSV

HTML

XLSX

PDF

Adjust color scale:

low

medium

high

default

0

0.05

0.3

Color Scheme:

RBG

RWB

YBB

Table Output Help

Mixture	apCAFs	B cells	C1QC+ Macrophages	CD8+ T cells	cDC1	cDC2	eCAFs	Endothelial cells	iCAF+
TCGA-KU-A6H7-06A-21R-A31N-07	0.02	0.00	0.00	0.05	0.00	0.02	0.00	0.00	0.00
TCGA-UF-A71A-06A-11R-A39I-07	0.38	0.00	0.00	0.00	0.00	0.00	0.07	0.03	0.00
TCGA-4P-AA8J-01A-11R-A39I-07	0.37	0.00	0.02	0.00	0.00	0.01	0.06	0.00	0.00
TCGA-BA-4074-01A-01R-1436-07	0.00	0.00	0.04	0.01	0.00	0.05	0.32	0.02	0.00
TCGA-BA-4075-01A-01R-1436-07	0.00	0.00	0.03	0.00	0.00	0.02	0.13	0.06	0.00
TCGA-BA-4076-01A-01R-1436-07	0.05	0.00	0.02	0.00	0.00	0.02	0.00	0.00	0.00
TCGA-BA-4077-01B-01R-1436-07	0.05	0.00	0.06	0.02	0.00	0.07	0.01	0.01	0.00
TCGA-BA-4078-01A-01R-1436-07	0.57	0.00	0.09	0.00	0.00	0.00	0.16	0.01	0.00
TCGA-BA-5149-01A-01R-1514-07	0.19	0.00	0.01	0.00	0.00	0.01	0.02	0.04	0.00
TCGA-BA-5151-01A-01R-1436-07	0.05	0.00	0.03	0.00	0.00	0.02	0.00	0.00	0.00
TCGA-BA-5152-01A-02R-1873-07	0.02	0.00	0.12	0.01	0.00	0.00	0.04	0.01	0.00
TCGA-BA-5153-01A-01R-1436-07	0.08	0.06	0.06	0.09	0.00	0.00	0.01	0.06	0.00
TCGA-BA-5555-01A-01R-1514-07	0.03	0.00	0.08	0.00	0.00	0.00	0.01	0.01	0.00
TCGA-BA-5556-01A-01R-1514-07	0.06	0.01	0.14	0.04	0.00	0.00	0.00	0.01	0.00
TCGA-BA-5557-01A-01R-1514-07	0.29	0.01	0.00	0.00	0.00	0.06	0.00	0.00	0.00